Simplify your complex SBRT fiducial marker placements!

CIVCO Cartridge Kits reduce the number of needle passes and minimize the chance for Edema or Pneumothorax



Insert and guide needle to desired location and remove the stylet.



Use CIVCO's pre-loaded cartridge to place fiducial marker into needle hub.



Return stylet to deploy marker. To create desired reference, reposition the needle to add additional markers.

Specialty Marker Kit Options*





PointCoil™ SBRT Kit

Trocar placement needle Qty 1 19 GA x 15 cm

Pre-loaded PointCoil cartridge Qty 3 0.8 x 5 mm marker

Additional Cartridge

Pre-loaded PointCoil cartridge Qty 1 0.8 x 5 mm marker


MTCTXPC0819153

MTCTXPC08185

MTCTXFC001





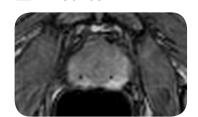

FusionCoil™ SBRT Kit

Pre-loaded FusionCoil needle Qty 1 18 GA x 20 cm; 1 x 5 mm marker

Pre-loaded FusionCoil cartridges Qty 2 1 x 5 mm marker

Additional Cartridge

Pre-loaded FusionCoil cartridge **Qty 1** 1 x 5 mm marker

MTCTXFC105185

*Not CE Marked

Global Sales Office Orange City, Iowa, USA 800.842.8688 I +1 712.737.8688

info@CivcoRT.com www.CivcoRT.com

