


Hip + Pelvic Positioning Devices

Carbon Fiber Bellyboard

The Carbon Fiber Bellyboard assists in reducing the volume of the small bowel that is irradiated in patients undergoing treatment in the pelvic region. The contoured opening allows for displacement of the small bowel, as the central support structure assists in compressing the small bowel forward into the opening.

- Carbon Fiber model features embedded fiducial crosshair to assist as a secondary setup reference
- Lightweight, rigid and indexable
- Lower part of board is adapted to support the patient's upper legs for maximum comfort
- Soft polyethylene padding ensures a comfortable position for the head

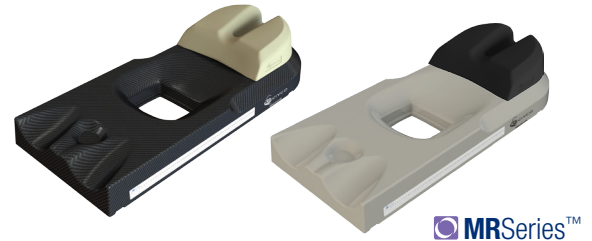

Carbon Fiber Bellyboard	125012
Prone Thorax Support (Cream) <i>sold separately</i>	301076
MR Bellyboard	126000
Prone Thorax Support (Black) <i>sold separately</i>	301102

Combifix™ 3

The Combifix 3 is a baseplate system providing comfortable positioning with enhanced immobilization when treating the pelvic region and other lower extremities. The system combines the Kneefix™ 3 and Feetfix™ 3 cushions with optional Elevation Blocks into a fixed and indexable position.

- Kneefix 3 height is adjustable in 5 cm increments with the use of Elevation Blocks
- Feetfix 3 is adjustable longitudinally in 8 indexable positions and 4 indexable angles
- Three-Pin Lok-Bar™ indexable and MR safe
- Feetfix Couch Support can be purchased separately
- Optional Shoulder Retractor can be attached in 5 positions to Feetfix Couch Support



Combifix 3 (includes Combifix Baseplate, Kneefix 3, Elevation Block 3, Feetfix 3 and Feetfix 3 Couch Support)	151004MR
Combifix 3 (includes Combifix Baseplate, Low Kneefix 2, Elevation Block 3, Feetfix 3 and Feetfix 3 Couch Support)	151005MR
Shoulder Retractors (MR Unsafe)	302072

Global Sales Office

Orange City, Iowa, USA
800.842.8688 | +1 712.737.8688

info@CivcoRT.com
www.CivcoRT.com

COPYRIGHT © 2018. CIVCO IS A REGISTERED TRADEMARK OF CIVCO MEDICAL SOLUTIONS. COMBIFIX, KNEEFIX, AND FEETFIX ARE TRADEMARKS OF CIVCO. ALL OTHER TRADEMARKS ARE PROPERTY OF THEIR RESPECTIVE OWNERS. ALL PRODUCTS MAY NOT BE LICENSED IN ACCORDANCE WITH CANADIAN LAW. 2018AP1348 REV. A

CIVCO[®]
Radiotherapy
a better way