

ThermFix™ Immobilization System

for Breast, Thorax, Abdomen, Hip, & Pelvic Treatment

The ThermFix Immobilization System uses a single sheet of thermoplastic that forms over the breast, thorax, abdomen, hip, or pelvis to aid patient positioning and reproducibility. The ThermFix Bars easily attach and index to the TotalRT™ All in One Platform or Universal Couchtop™ and provide 1 cm indexing.

Features & Benefits

- Bar design provides a larger and more homogeneous treatment area compared to a baseplate
- Multiple length and width indexing options accommodate many patient sizes
- Allows supine or prone positioning
- The inferior bar can be removed during patient load/unload to provide safer transitions
- Vac-Lok™, RediFoam™ Cushions, and Wing Boards™ can be integrated for added comfort and support
- MR Safe, CT, and treatment compatible



ThermFix Base Package

Includes: Superior Bar, Inferior Bar, Spacing Tool (for breast, thorax, abdomen, hip and pelvic treatment with no central docking)

THRFX100



ThermFix Pelvic Package

Includes: Superior Bar, Hip & Pelvic Bar, Spacing Tool (for hip and pelvic treatment with central docking)


THRFX101



ThermFix Complete Package

Superior Bar, Inferior Bar, Hip & Pelvic Bar, Spacing Tool

THRFX102

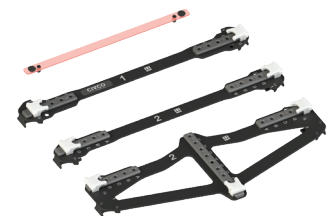


ThermFix™ Thermoplastics

ThermFix Thermoplastics are available with reloadable or disposable handles with solid, standard perforated, or IMRT Reinforced™ perforation styles. The handles and groin piece secure to the ThermFix Bars, providing an indexable setup.


Part Number	Perforation	Handle Type	Thickness
MTHFX1822S	Solid	Reloadable	3.2 mm
MTHFX18222.4S	Solid	Reloadable	2.4 mm
MTPLHA0201	Solid	Disposable*	3.2 mm
MTHFX1822	Perforated	Reloadable	3.2 mm
MTHFX18222.4	Perforated	Reloadable	2.4 mm
MTPLHA0202	Perforated	Disposable*	3.2 mm
MTHFI1822	IMRT Reinforced	Reloadable	3.2 mm
MTHFI18222.4	IMRT Reinforced	Reloadable	2.4 mm
MTPLHAI0203	IMRT Reinforced	Disposable*	3.2 mm

MTPLHA01	Reloadable Handles (x2) & Groin Piece
MTPLVCSUB8	Reusable Groin Piece

*Groin piece sold separately

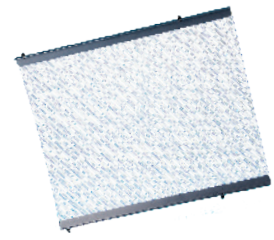

Standard Perforated Thermoplastic

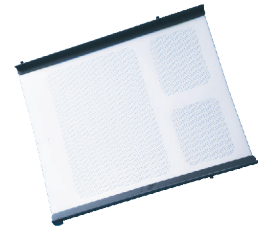

IMRT Reinforced Thermoplastic



Global Sales Office

Orange City, Iowa, USA
800.842.8688 | +1 712.737.8688

info@CivcoRT.com
www.CivcoRT.com

